MADE IN CHINA

RECOVERYM8

USER MANUAL

Recoverym8 Smartcup

SMART CUPPING

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ENGLISH MANUAL

SMART CUPPING

The Smart cupping is a dynamic auto cupping therapy device with a revolutionary built-in automatic air compression system.

It is a game-changer in cupping therapy with its dynamic cupping technology. There are five levels of custom intensity with five modes to adjust for different massage areas and levels of activity in each person. It also integrates red light and is effective in relieving pain and tension in the muscle. It is an excellent theracup addition for massage use and has a compact design, making it convenient for home/office/travel use.

It is also a scraping massage tool with innovative adjustable suction power that helps relieve tension in the muscles when used with some lotions. It also helps improve blood circulation and lymphatic drainage to reduce bloating.

50
dB

Quieter motor

0.2
KG

Light weight

5
W

Power

3
H

Battery life

MAIN HIGHLIGHT FEATURES

Dynamic Cupping

Smart cupping exceeds traditional cupping therapy sets with its dynamic suction. The alternation of suction and release to simulate authentic massage boosts efficient blood circulation while speeding up the recovery process.

Red Light Therapy

The smart cupping has added Red Light Therapy (RLT) during the cupping therapy, which helps relieve pain and inflammation to achieve a more effective treatment.

Exposing body parts under red light also yields faster results in healing the skin and muscle tissues and helps you relax better.

Personalized Massage

Faster or slower? Stronger or less suction? Not a problem, Recoverym8 Smartcup allows you to set your favorite modes and comfortable suction power. Make the cupping therapy more amusing.


Guasha / Scraping

Suction mode number 5 is the scraping massage mode, excellent for people who want to relieve tension in the muscle, boost circulation, and improve lymphatic drainage.

Weight: **0.44 lbs / 200g**
Dimensions: **3.4x3.4x2.83 inches**
Charging voltage: **5V-2A**
Suction intensity: **-60~0 kPa**

Charging time: **2.5 Hrs**
Wavelength: **630-670nm**
Power rating: **5W**
Battery capacity: **3.7V / 1800mAh**

- 1 Suction Mode Indicator
- 2 Suction Mode Control
- 3 Suction Movement Light
- 4 Suction Power Control
- 5 Suction Power Indicator
- 6 Power Button
- 7 Charging Indicator
- 8 Charging Port
- 9 Lanyard Hole
- 10 Hard Plastic Cup
- 11 Air Inlet
- 12 Red Light Zone

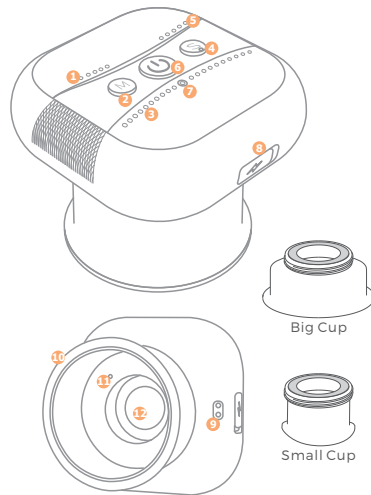QUICK START GUIDE

Before Use

Please read precautionary statements first. Do not use the Recoverym8 Smartcup if not suitable.

Moisturize your skin with lotions or massage oils for a more comfortable cupping experience.

Please clean the product with a wet rag before and after use.

Tighten the cup on the main unit and check if the machine is complete.

Step 1: Power Button

Press and hold to power on/off the device.

(Suction is on by default once power is on)

Press the power button to pause/continue.

Auto turn-off by default after 12-min continuous working.

Step 2: Set Suction Mode

Mode 1 is on by default once power is on.

Press the "M" button to switch between 5 different suction modes.

(Indicated by different LED green lights)

Mode 1 - Massage Mode Mode 2 - Radical Mode Mode 3 - Relaxing Mode

Mode 4 - Cupping Mode Mode 5 - Scraping Mode

Blue LED lights will flash to simulate the pumping dynamic motion.

Step 3: Set Suction Power Level

Level 1 intensity is on by default once power is on.

Press "S" button to the custom suction intensity from level 1 to level 5.

The indicator red light shows the current intensity level accordingly.

FRIENDLY TIP

There's a tactile dot on the "S" button. Use it as a blind spot indicator when it is on your back so that you can find it quickly to set the suction power level. On the contrary, the "M" button without a blind spot indicator is the suction mode setting.

ACCESSORIES USE

Filter Cotton Ring

Filter cotton ring is designed to block out most granular or liquid impurities for internal cleanliness and maintenance.


Please check to replace it based on the frequency of use.

(Presented three more in the standard kit)

Replacement steps


Step 1: Twist the plastic cup off.

Step 2: Replace the cotton ring with a new one and then tighten the cup.

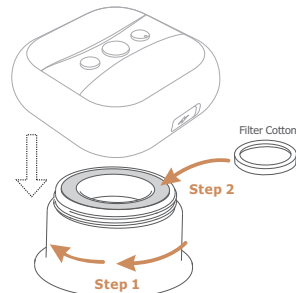
Lanyard

Please put the lanyard on and hold it to prevent the product from falling and getting damaged.

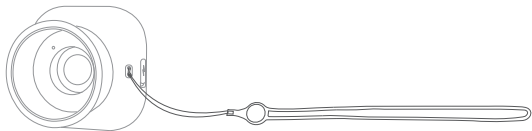
Type-C charging cord

Standard type-C charging cord is provided.

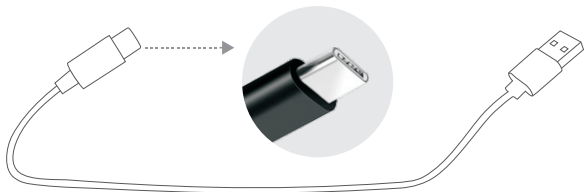CHARGING GUIDE

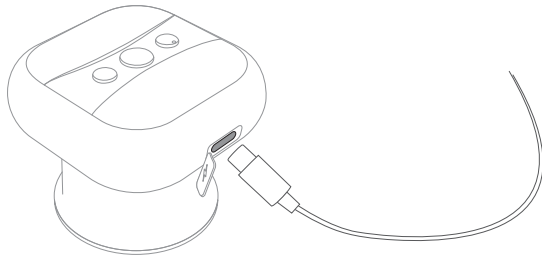1. To charge the Smart cupping, please connect USB A end of the charging cord to any adapter that supports voltage 5V, current from 1A~2.5A, and the other end USB C to the product.
2. The charging time is about 2.5 hours. The charging indicator light is red while charging and turns blue when fully charged. Due to environmental conditions and different charging adapters, the actual charge time may vary.
3. Please Keep the charging cord dry and clean, and do not overcharge for more than 8 hours.
4. Please do not use the device when it is charging.
5. When the battery voltage is too low, it automatically shuts down after beeping ten times.


CLEANING & MAINTENANCE

Cleaning

Please turn off the machine before cleaning and use a dry non-woven cloth to clean the machine body.

Do not clean the machine body with acid and alkaline cleaners, alcohol, thinner, gasoline, nail polish remover, and other solvents, stiff brush, or stainless steel sponge.

Do not use liquids such as water to clean the machine body directly.

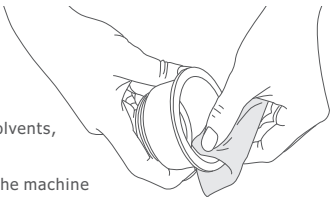Do not use and store this product in the following places

Store this product in a safe and dry place.

Do not leave the product idle, with moisture remaining.

Do not use or store this product in places susceptible to vibration, under direct sunlight, with high temperature and humidity, with water, near strong electromagnetic fields or anything explosive, etc.

SAFETY INSTRUCTIONS

Danger

- Do not use it with medical implants (e.g. cardiac pacemakers), life-sustaining medical electronic devices (e.g. artificial heart and lungs), medical electronic equipment (e.g. electrocardiograph).
- Keep out of reach of children.

Warnings

- Do not apply to wounds, scars or areas near the heart.
- If discomfort occurs, stop using it and inform your healthcare provider immediately.
- Attach the plastic cup evenly to the skin.
- Recommend use on one person only. To apply it to another person, please first turn off the machine and clean the plastic cup.
- Avoid using it on the same body part for too long (device shuts off at 12 minutes) to prevent irritation and possible skin damage.
- Do not use over an area where the skin is bruised from previous use.
- Do not compress internal organs of the body with the device to prevent accidents or physical discomfort.
- Do not use a chipped or damaged plastic cup.
- Do not use it for medical effects.
- This product contains lithium batteries. Please do not disassemble, modify, immerse it in water or throw into fire to prevent damage and explosion.

- Please obtain approval from the doctor if any of the following apply:
allergies; infectious diseases; body temperature above 38°C; during the medication period; allergic to metals and red light; disorders of the heart, cranial nerves, blood pressure; disorder of perceptions or the skin; sensory disorders due to other high-level peripheral circulatory disturbances of diabetes or other diseases; cannot express one's own will freely; being treated by a doctor or feel unwell.
- It's prohibited to use this product if any of the following apply:
 1. Women during pregnancy or menstrual cycle;
 2. Patients with malignant tumours or local damage and ulceration, traumatic fractures, varicose veins, acute suppurative inflammation, and loss of skin elasticity;
 3. Severe heart disease, heart failure, respiratory failure and severe edema;
 4. Patients with poor coagulation mechanism, spontaneous bleeding or bleeding after injury;
 5. People with severe skin allergies or skin infectious diseases;
 6. Severe nervousness, high fever, convulsions, spasms, manic restlessness;
 7. Area around heart, sensory organs, private areas or large blood vessels on the body.

Caution

- Please use this product following the instructions.
- It's normal for the product to run hot after using for a long time.

- Redness or bruise might occur after using the product, which is normal and usually disappears in 3- 5 days (depending on each person's skin type).
- Do not plug or unplug the USB charging cable with wet hands to prevent electric shock or accidents.
- Please do not use it on children or people without expressive ability to prevent accidents.
- Do not use this product while driving.
- Do not stick the product on the affected area for a long time to prevent skin inflammation.
- The battery will run off if idle for a long time. Please follow the instructions and recharge before use.
- The product contains batteries. Please recycle them when the service life is over.
- Best use by people over 16 years old.

WARRANTY STATEMENT

This product enjoys a two-year warranty. The starting date of the warranty period is based on the product invoice date. If the product fails due to defects in materials or design within two years, the company will repair or replace new ones for free, except in the following cases:

1. Improper use or damage to the product by the user.
2. Unauthorized disassembly and maintenance of the product.
3. Damage due to using gifts or accessories that are not officially offered.
4. The model included in the invoice does not match the model of the repaired product or has been altered.
5. Gifts and accessories.
6. Out of warranty period or lack of valid invoice.
7. Damage caused by force majeure.

FAQ

QUESTIONS	ANSWERS
Can not turn on the device?	1. Hold the power button over 1.5 seconds. 2. Charge the machine for 2 hours and try again.
Weak suction	1. Re-install and tighten the plastic cup. 2. Check if the air inlet/outlet is blocked. 3. Check if the power button cannot be reset. 4. Check if the filter cotton needs replacement.
Can not charge?	1. Try again with a different socket. 2. Pull out the adapter and try again to make sure it is in place. 3. Check and clean the adapter and charging port if there is any dirt or debris. 4. Replace the adapter or charging cable.
Why blue light when power is on?	It turns to red light when the machine suction is applied to your body.
Can I get more different size cups?	All 3 different size cups are presented in the standard kit.
What if I lose my charging cord?	Smart cupping uses a regular charging cord with USB A to USB C, so that the charging cord can be replaced with any such kind if you lose the original one. Or buy a new one from us.
How long should I cup per spot?	It's suggested to cup from 3 to 12 minutes per spot. You could also use a timer to keep track of the time.
How do I turn it off while it's running?	Just push the on/off button to stop it.
How long does it take for the cupping marks to fade away?	According to our study, it normally takes around 3-5 days for the cupping marks to fade away. It might take a bit longer for some individuals.

Do not disassemble the machine and please contact us if the issue can not be solved.

Accessories

Lanyard * 1

Charging Cable * 1

Filter Cotton * 3

Big Cup * 1

Small Cup * 1